

MERCATOR nitrile blue

RP30007001-05

The instruction below should be used in conjunction with detailed information on the packaging.

Short description of the product

Nitrile protective gloves, powder-free, non-sterile for disposable use, for minimal risk

Full description of the product

Reference number: : RP30007001-05
Raw material : nitrile
Cuff : beaded
Colour : blue
Shape : ambidextrous, fitting either hand
Size range : XS (5-6), S (6-7), M (7-8), L (8-9), XL (9-10)
Quantity in packaging : 100 pcs. by weight
Shelf life : 3 years (from the date of manufacturing)

Storage instructions


It is recommended to store the gloves in dry place, in the temperature of 5-35°C and to protect them against direct sunlight and fluorescent light.

It is recommended to protect gloves against humidity.

Keep the gloves in a distance of not less than 1m from heating devices, sources of fire and ozone.

Do not keep in direct vicinity of solvents, oils, fuels and lubricants.

Food contact

Gloves are marked with food contact symbol  and comply with the requirements of Regulation (EU) No 10/2011, European Regulation (EC) No 1935/2004 and with Regulation (EC) No 2023/2006 on Good Manufacturing Practice.

PPE classification & compliance

Gloves are classified as Personal Protective Equipment – category I for minimal risk according to the Regulation 2016/425.

Compliance with standards: EN 420:2003+A1:2009.

Declaration of Conformity can be found under below web address: <https://en.mercatormedical.eu/>



Intended use

These are non-sterile protective gloves, intended for use in hygienic and sanitary works in public places, for maintaining cleanness in facilities, utility rooms and home, professional cleaning services using cleaning materials of weak action, gardening, agriculture, food processing and food contact, for use in household and mechanical works with risk only superficial mechanical injury. Classified as Personal Protective Equipment category I for minimal risk. Their design and labeling corresponds to the requirements of the Regulation 2016/425 on Personal Protective Equipment. Gloves should be used solely according to their intended purpose.

Precautions and indications for use

Dry hands before putting the gloves on. Before usage, inspect the gloves for any defect or imperfections. If the gloves get punctured, torn or broken during their use, take them off and put on the new ones. Avoid using gloves dirty in the inside as they may cause irritation leading to skin inflammation or more serious damages. The gloves should not be used in contact with open fire and to protect against any sharp tools. The gloves are not intended for welding, electric shock protection, ionizing radiation or from the effect of hot or cold objects.

Components / hazardous components

Some gloves may contain components known to be a possible cause of allergy for person allergic to them, who may develop contact irritation and/or allergic reaction. In the case of an allergic reaction, seek medical assistance immediately.

Disposal

Used gloves should be disposed of in accordance with local regulation. Gloves should be buried or burned under controlled conditions.


Manufacturer

MERCATOR MEDICAL S.A.
ul. H. Modrzejewskiej 30
31-327 Cracow, Poland
www.mercatormedical.eu

Symbols used on the packaging

	Medical device		Personal Protective Equipment		Powdered gloves
	Do not re-use / gloves are intended for single use		Keep away from moisture, store in a dry place		Powder free gloves
	Non-sterile gloves		Keep away from solar and fluorescent light		Presence of polymer coating on the inner surface of the glove
	Lot / batch number		Temperature limitation / gloves store in temperature 5-35°C		Presence of cosmetic coating on the inner surface of the glove
	Catalogue number		Keep away from ozone		Presence of external texture on the glove
	EU Authorised Representative, symbol should be accompanied by name and address of Authorised Representative		Date of manufacture		Gloves made from nitrile
	Expiry date		Manufacturer, symbol should be accompanied by name and address of Manufacturer		Gloves made from vinyl
	Marking of gloves protecting against bacteria and fungi		Food contact symbol (article is suitable for food contact, for details check the instruction for use)		Raw material – natural rubber latex
	Marking of gloves protecting against viruses, bacteria and fungi		Package made from paper, qualify for recycling		50 gloves by weight
	Marking of type A chemical resistant gloves. Six tested chemicals shall be identified by their code letter under pictogram		Package is treated as municipal waste		100 gloves by weight
	Marking of type B chemical resistant gloves. Three tested chemicals shall be identified by their code letter under pictogram		Consult instructions for use		200 gloves by weight
	Marking of type C chemical resistant gloves. One tested chemicals shall be identified by their code letter under pictogram		Information in other languages inside		Do not use, if package is damaged
	Indicates compliance with the requirements of Russian market		Indicates compliance with the requirements of Ukrainian market		

HOW TO PUT THE GLOVES ON?

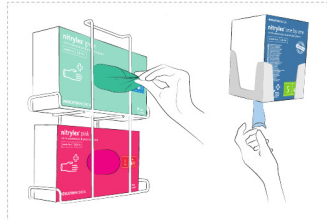

Take one glove from the box pulling by the cuff.



Hold the glove by the cuff with the same hand you used to pull it from the box and put the other hand into the glove without touching the working surface.



Stretch the glove on the right hand.



Take another glove from the box with the hand wearing a glove.




Put the glove on the hand and smooth it out.




Done!


HOW TO TAKE THE GLOVES OFF?

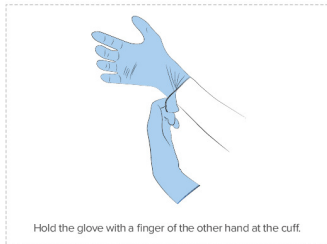


Hold the glove with a finger of the other hand at the cuff.

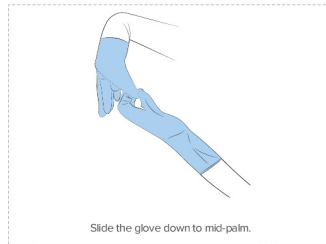


Slide the glove down to mid-palm.

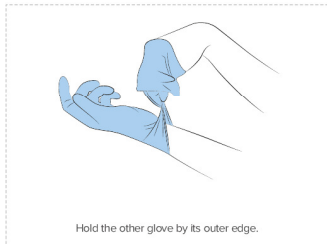

Hold the other glove by its outer edge.

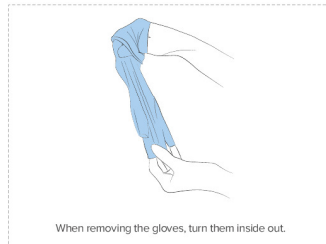

When removing the gloves, turn them inside out.

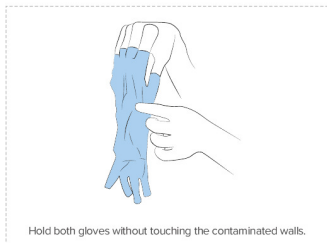

Hold both gloves without touching the contaminated walls.



Removed gloves should be treated as healthcare waste.